

# Meditherapy®

Accessories  
For Endoscopy




ENDOBITE COATED BIOPSY FORCEPS

# ENDO Bite

## SINGLE-USE BIOPSY FORCEPS

DESIGNED AND PRODUCED BY  
**medwork**®

By utilising the same pantograph hinge mechanism as reusable biopsy forceps to ensure optimum opening of the cups and optimal bite force, the EndoBite biopsy forceps guarantee an excellent quality of biopsy even when the endoscope is tightly angulated.

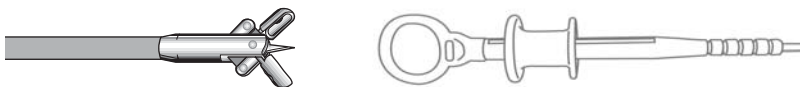

### FEATURES

- Pantograph opening of the cups with high precision mechanism, traditionally used for reusable forceps, provides added bite force for increased biopsy quality
- Sharp edged cups to enable 'cut' instead of 'tear' biopsy
- 'Channel Guard' low friction sheath is a unique fluoropolymer coating offering a very low coefficient of friction, facilitating easy passage of the biopsy forceps within the instrument channel. The lower friction sheath may also reduce wear to the biopsy channel wall
- Fenestrated cups for an increased biopsy volume
- Colour coded to differentiate between different forceps models
- Spiked forceps are delivered with a protective distal cap
- Traceability : 3 detachable stickers accompany each pair of forceps



SINGLE USE BIOPSY FORCEPS					Class IIa CE0123
Part number	Compatible channel	Length (cm)	Spike	Colour code	Pack size
BIO1-C4-23-180	2.8 – 4.2	180	No	Green slider	10
BIO1-C4-23-230	2.8 – 4.2	230	No	Red slider	10
BIO1-D4-23-180	2.8 – 4.2	180	Yes	Green slider/ Black insert	10
BIO1-D4-23-230	2.8 – 4.2	230	Yes	Red slider/ Black insert	10
BIO1-C4-18-120	2.0 – 3.2	120	No	Blue slider	10
BIO1-C4-18-180	2.0 – 3.2	180	No	Blue slider	10
BIO1-C4-18-260	2.0 – 4.2	260	No	Green slider	10