
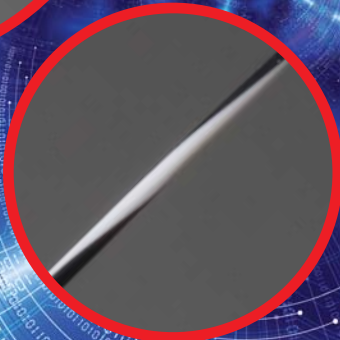


# Meditherapy®

Accessories  
For Endoscopy

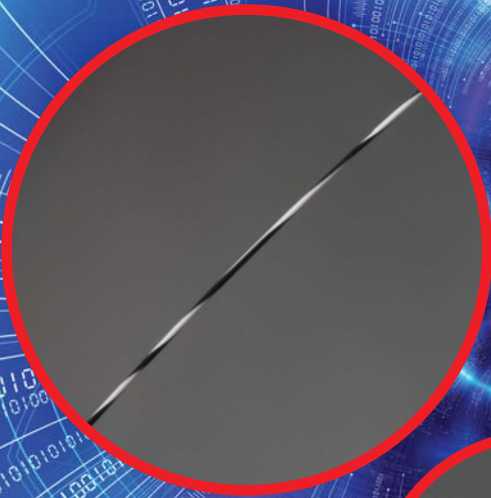

# WIR1 SERIES

## ERCP NEG GUIDEWIRES

DESIGNED AND PRODUCED BY



The WIR1 Series of NEG ERCP guidewires are especially designed to aid endoscopic visual movement of the wire whilst providing excellent handling characteristics to meet different requirements and preferences.



### FEATURES

- The tungsten impregnated 5cm floppy tip provides enhanced cannulation and greater fluoroscopic visibility
- The Nitinol core is kink resistant and thus aids device exchange
- The NEG wire incorporates 1cm black and white visual markers to help detect inadvertent wire movement
- The insulated PTFE jacket conforms to the ANSI/AAMI HF 18 standard enabling a sphincterotomy knife to cut over wire without the need to remove it
- The wire comes in standard, soft, J-Tip and 0.025" and 0.035" versions for all eventualities



STERILE

### SINGLE USE NEG GUIDEWIRES

Class IIa CE0123

Part number	Size (inches)	Length (cm)	Type	Pack size
WIR1-F5-25-400	0.025"	400cm	Standard	2
WIR1-F5-35-400	0.035"	400cm	Standard	2
WIR1-G5-35-400	0.035"	400cm	Extra soft	2
WIR1-H5-35-400	0.035"	400cm	J-Tip	2