


# Meditherapy<sup>®</sup>

Accessories  
For Endoscopy

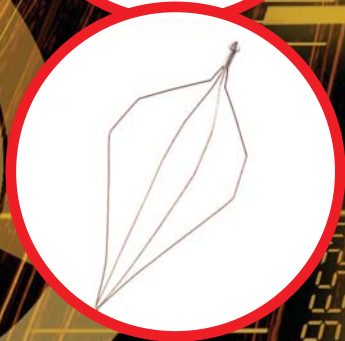

# BAS1 AND FOR1

## GRASPING FORCEPS AND RETRIEVAL DEVICES

DESIGNED AND PRODUCED BY  
**medwork**<sup>®</sup>

The ERCP BAS1 basket design provides a 'rounder' shape to facilitate opening and stone capture. In addition the rat tooth forceps offer exceptional grasping ability for biliary stent retrieval.

The three prongs of the tripod forceps have atraumatic ring ends to prevent polyp damage and an opening width of 25mm enables the capture of a wide variety of different sized polyps.

### FEATURES

- The superior 16 sided basket composed of 4 multifilament wires with a diameter of 0.4 mm enables excellent opening dynamic as well as optimal rounder shape retention. These features facilitate easier stone capture
- Baskets incorporate an injection port for contrast injection
- All baskets have a built in safety point at the proximal end enabling emergency lithotripsy to take place if required
- Unique 'Twister' basket rotates 130° upon closure facilitating stone capture
- Rat tooth forceps with wide jaw opening are purpose designed for biliary stent removal
- Tripod polyp retrieval forceps incorporate rounded prong ends for atraumatic polyp retrieval. The thin flexible prongs prevent crush artefact while the short throw handle provides perfect ergonomic control



STERILE

### SINGLE USE EXTRACTION BASKETS

Class IIa CE0123

Part number	Channel size	Length (cm)	Diameter of basket (mm)	Type of Basket	Pack size
BAS1-A2-30-23-200	2.8mm or larger	200	30	4 wire hexagonal	2
BAS1-H1-20-27-200-E35	2.8mm or larger	200	30	Wire guided 4 wire hexagonal	1
BAS1-K2-25-23-200	2.8mm or larger	200	25	4 wire spherical twister basket	2

STERILE

### REUSABLE RAT TOOTH RETRIEVAL FORCEPS

Class IIa CE0123

Part number	Channel size	Length (cm)	Type of grasp	Opening width	Pack size
FOR2-A1-23-230	2.8mm or larger	230	Rat tooth	12mm	1



STERILE

### SINGLE USE TRIPOD RETRIEVAL FORCEPS

Class IIa CE0123

Part number	Channel size	Length (cm)	Type of grasp	Opening width	Pack size
FOR1-H3-23-230	2.8mm or larger	230	Tripod type	25mm	5