


Meditherapy®

Accessories
For Endoscopy

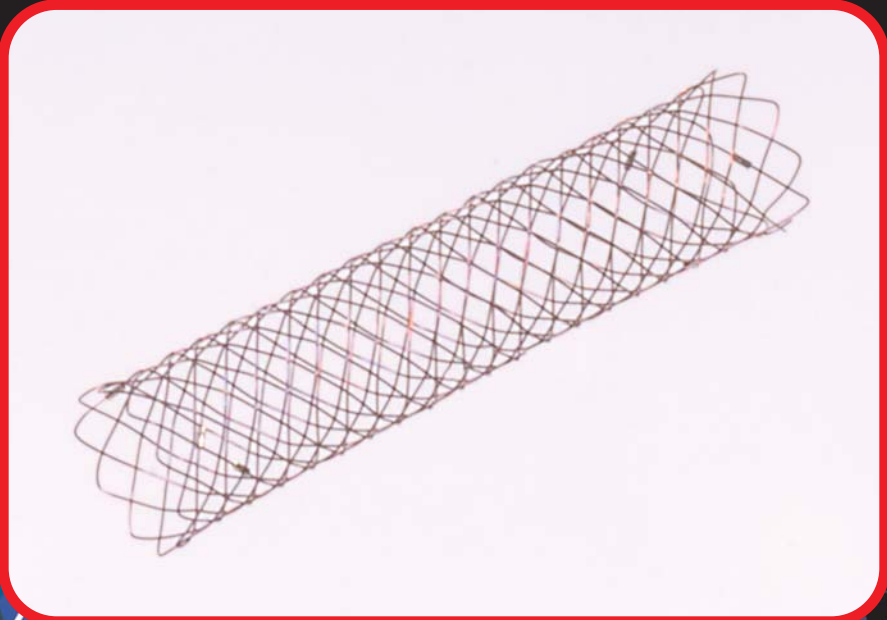

IMP BILARY STENT RANGE

IMP1 SERIES

BILIARY METAL SELF-EXPANDING STENTS

DESIGNED AND PRODUCED BY



The outstanding IMP nitinol metal stent range sets new standards in stent technology. The flexible, self-expanding stents have exceptional radiological visualisation characteristics, aiding in precise positioning of the stents.

FEATURES

- Nitinol construction ensures that the maximum diameter and radial force are reached at body temperature
- The superb expansion capabilities provides excellent opening ability even through difficult or torturous strictures
- Atraumatic edges are rounded to minimise possible mucosal damage
- IMP stents have several platinum–iridium alloy markers in the middle and at both ends of the stent to aid accurate placement
- IMP stents have exceptional flexibility to adapt to different anatomical conditions
- Uncovered and partially covered stents coupled with a range of diameters and lengths ensure all clinical situations are catered for
- The slim 7 and 8fr Insertion system diameter aids accurate placement of the stent.
- IMP stents are 0.035" guidewire compatible
- Two radiopaque external makers are supplied with each stent



STERILE

IMP1 BILIARY EXPANDABLE STENTS		Class IIa CE0459
Part number	Description	Pack size
IMP1-B2-10-040	Uncoated. Diameter 10mm, Length 40mm. 0.035" guidewire compatible	1
IMP1-B2-10-060	Uncoated. Diameter 10mm, Length 60mm. 0.035" guidewire compatible	1
IMP1-B2-10-080	Uncoated. Diameter 10mm, Length 80mm. 0.035" guidewire compatible	1
IMP1-B3-10-040	Covered. Diameter 10mm, Length 80mm. 0.035" guidewire compatible	1
IMP1-B3-10-060	Covered. Diameter 10mm, Length 60mm. 0.035" guidewire compatible	1
IMP1-B3-10-070	Covered. Diameter 10mm, Length 70mm. 0.035" guidewire compatible	1