
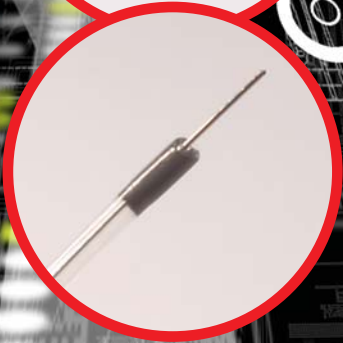


# Meditherapy®

Accessories  
For Endoscopy

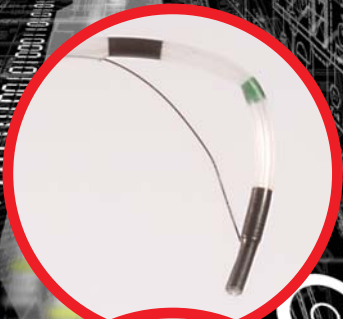

# PAP1 SERIES

## SPHINCTEROTOMY KNIVES

DESIGNED AND PRODUCED BY  
**medw<sup>ork</sup>**®

The precise manufacturing process of the PAP1 range ensures cutting takes place at the 11 o'clock position. These advances in sphincterotome technology will aid consistent and repeatable cannulation to ensure safer sphincterotomy.



### FEATURES

- Precise wire orientation always at the 11 o'clock position provides extra safety during sphincterotomy
- 4.5fr smooth tapered catheter aids cannulation even for the triple lumen knife
- 0.035" guide wire compatible for ease of use
- Coloured cutting wire markings provide visual identification of the precise length of cutting wire
- Three ring ergonomic handle design provides an excellent tactile feel during knife bowing
- Pre-curved for easier cannulation
- Tantalum radiopaque metal marker for clear identification of tip placement
- Connectable to majority of monopolar electro-surgical units



STERILE

SINGLE USE SPHINCTEROTOMY KNIVES					Class IIa CE0123	
Part number	Guide wire compatibility	Nose Length(mm)	Wire Length (mm)	Type	Pack size	
PAP1-C8-20-35-OL	0.035"	6	20	Triple lumen knife	1	
PAP1-C8-30-35-OL	0.035"	6	30	Triple lumen knife	1	
PAP1-K1-06-OL	-	6	-	Single lumen needle knife	1	
PAP1-K2-06-35-OL	0.035"	6	-	Double lumen needle knife	1	
PAP1-M2-20-35-OL	0.035"	6	30	Double lumen BII knife	1	