

Meditherapy®

Accessories
For Endoscopy




POL1 SERIES

POLYPECTOMY SNARES

DESIGNED AND PRODUCED BY

medwork®

The reduced throw handle of the POL1 range of polypectomy snares together with the 3-ring design, enhances tactile response while making the opening of the snare feel more comfortable for the nurse.

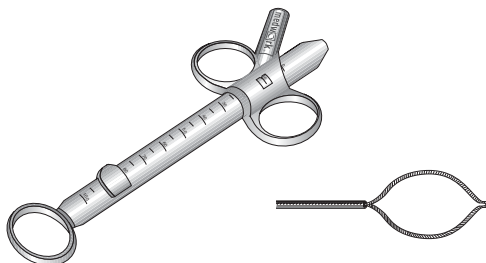

FEATURES

- Precisely engineered three ringed handle for added user comfort. Includes an integral handle gauge which permits the user to estimate the size of the captured polyp or pedicle and help ensure that normal tissue has not also been trapped within the loop
- Multifilament small diameter braided wire for precise polypectomy and optimal coagulation
- Rotatable version for difficult to capture polyps
- Unique 'Combi' snare combines both a large and small loop in one device. Assists in inventory management and ensures the correct size of snare is used every time
- Unique 'Resection Master' flat ribbon snare for EMR and sessile polyps
- Teflon™ sheath for optimal protection of the endoscope, resistance to high temperatures and excellent insertion through the instrument channel
- Inner centering tube, to help reduce the throw of the snare
- Connectable to majority of monopolar electrosurgical units



STERILE

SINGLE USE POLYPECTOMY SNARE

Class IIa CE0123

Part number	Compatible channel (mm)	Length (cm)	Loop width (mm)	Loop shape	Pack size
POL1-B7-30-23-220-OL	2.8 – 4.2	220	30	Oval	10
POL1-B7-20-23-220-OL	2.8 – 4.2	220	20	Mini Oval	10
POL1-B8-30-23-220-OL	2.8 – 4.2	220	30	Rotatable Oval	5
POL1-C5-30-23-220-OL	2.8 – 4.2	220	30	Flat Oval	5
POL1-D1-30-23-220-OL	2.8 – 4.2	220	15 and 30	Combi Oval	10
POL1-F1-30-23-220-OL	2.8 – 4.2	220	30	Hexagonal	10