

# Meditherapy®


Accessories  
For Endoscopy



# UNO-BRUSH

## CYTOLOGY BRUSH

Excellent manufacturing quality, and innovative design features make the Uno-Brush one of the very best cytology brushes available.

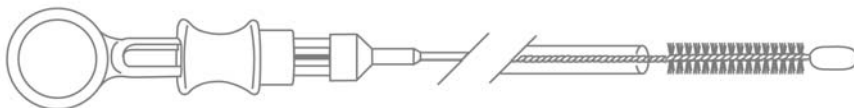

### FEATURES

- Occlusive distal tip protects brush from cell contamination during insertion and withdrawal through the scope
- Ergonomic handle is easily operated in one hand
- Rounded distal tip avoids damage to the instrument channel
- Includes two detachable stickers for traceability
- Single-use to reduce cross contamination

2

STERILE

### SINGLE USE CYTOLOGY BRUSH

Class IIa CE0459

Part Number	Description	Ø catheter (mm)	Length (cm)	Ø brush (mm)	Brush length (mm)	Pack size
AS0372512	Brush for bronchoscope	2.0 – 2.2	100	2.5	11	5
AS0373016	Brush for gastroscopy	2.8 – 4.2	160	3.0	11	5
AS0373022	Brush for colonoscopy	2.8 – 4.2	220	3.0	11	5