

# Meditherapy®

Accessories  
For Endoscopy




UNO-CLEAN CLEANING BRUSHES

# UNO-CLEAN

## CHANNEL CLEANING BRUSHES

The Uno-Clean single-use range of cleaning brushes are available in variety of designs for majority of endoscope channel sizes. Each brush has thick nylon bristles together with full rounded head and tapered distal end to enhance passage through scope. The range includes single ended, double ended and combination style products.

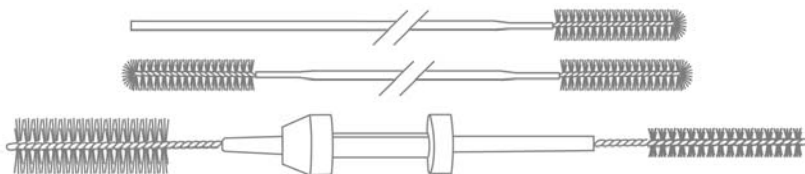

### FEATURES

- Brush diameter, length and flexibility ensure effective and thorough cleaning
- Full range available, including combination brush
- Brush head contains bristles around tip for frontal cleaning in angled areas of scope, without risk of damaging the channel
- Tapered distal end and reinforced metallic core, ease passage through angled channels and increases brush efficiency
- Traceability: 2 detachable stickers accompany each cleaning brush
- Separate Port/Valve cleaning brush available for cleaning inside channel inlets and outlets
- Colour coded for immediate identification of brush to be used



SINGLE USE CLEANING BRUSH							Class I CE	
Part Number	Number of brush heads	Colour	Compatible channel (mm)	Sheath length (cm)	Brush diam (mm)	Brush length (mm)	Pack Size	
AN1050110	1 for inst channel	White	1 - 1.7	100	2	8	50	
AN1050138	1 for air/water channel	White	1 - 1.7	350	2	8	50	
AN1053010	1 for inst channel	Green	2 - 2.2	100	3	18	50	
AN1053018	1 for inst channel	Green	2 - 2.2	180	3	18	50	
AN1055023	1 for inst channel	Yellow	2 - 4.2	230	5	20	50	
AN1065023	2 for inst channel	Yellow	2 - 4.2	230	3	20	50	
AN2085023	Combination brush, 1 for inst channel and integral port cleaning brush	Yellow	2 - 4.2	230	5&11	20	50	
KJ2085023	Combination brush, 1 for inst channel and includes integral port cleaning brush	Yellow	2 - 4.2	230	5&11	20	50	
AN1061110	Port and valve lumen brush	Blue		11	5&11	20&35	50	