


Meditherapy®


Accessories
For Endoscopy

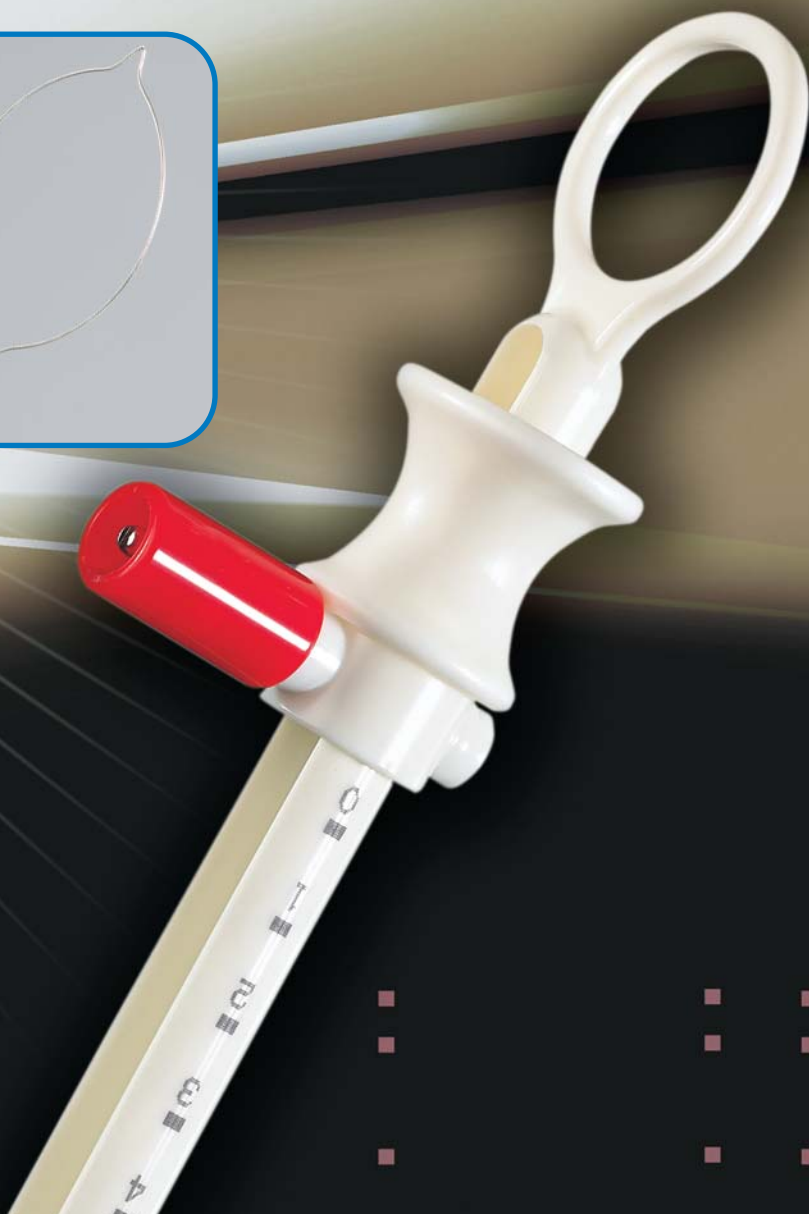

UNO-LOOP

POLYPECTOMY SNARE

The Uno-Loop single-use range of small diameter, multi-braided wire, snares provide precision polypectomy and effective coagulation. In addition, handle is fully IEC 601-2-18 compliant, providing total electrical safety for user.

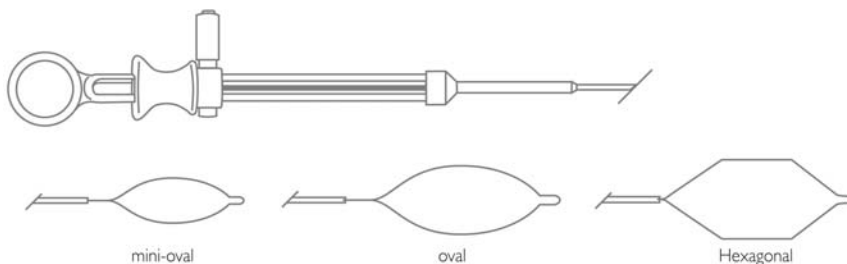

FEATURES

- Multifilament small diameter braided wire for precise polypectomy and optimal coagulation
- Fully conforms to IEC 601-2-18 electrical safety standard by providing insulation over complete length of snare (including entrance into the handle) to protect user
- Teflon™ sheath for optimal protection of endoscope, resistance to high temperatures and excellent insertion through the endoscope channel
- Inner centering tube, to reduce throw of snare
- Connectable to majority of monopolar electro-surgical units
- Colour-coded to enable immediate identification of loop shape



STERILE

SINGLE USE POLYPECTOMY SNARE						Class IIa	CE0459
Part Number	Compatible channel (mm)	Length(cm)	Loop width (mm)	Loop shape	Colour code	Pack size	
AS1525023	2.8 – 4.2	230	30 x 50	Oval	Blue	5	
AS1532523	2.8 – 4.2	230	15 x 25	Mini-oval	White	5	
AS1545023	2.8 – 4.2	230	30 x 50	Hexagonal	Grey	5	